

PATENT ABSTRACTS OF JAPAN

(11)Publication number : 2000-185023

(43)Date of publication of application : 04.07.2000

(51)Int.Cl.

A61B 5/055
G01R 33/48

(21)Application number : 10-363988

(71)Applicant : GE YOKOGAWA MEDICAL
SYSTEMS LTD

(22)Date of filing : 22.12.1998


(72)Inventor : GOTO TAKAO

(54) MR IMAGING METHOD, PHASE SHIFT MEASURING METHOD, AND MR IMAGING DEVICE

(57)Abstract:

PROBLEM TO BE SOLVED: To prevent the influence of fluctuation of residual magnetization resulting from changes in gradient magnetic field pulses.

SOLUTION: In an FSE sequence, reset pulses g_{yslj} and g_{ysrj} which have amplitude equal to maximum amplitude used on a phase axis are inserted on the phase axis in front of and behind an inverted pulse P_j . Therefore, residual magnetization prior to application of phase encoding pulses is always constant irrespective of the history of application of other gradient magnetic field pulses. Thus, the influence of fluctuation of residual magnetization resulting from changes in gradient magnetic field pulses can be prevented.

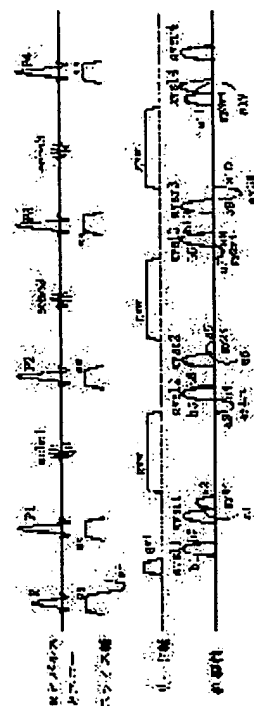

LEGAL STATUS

[Date of request for examination]

23.05.2000

[Date of sending the examiner's decision of rejection]

[Kind of final disposal of application other than the examiner's decision of rejection or application converted registration]

[Date of final disposal for application]

[Patent number]

3382868

[Date of registration]

20.12.2002

[Number of appeal against examiner's decision of

FB

rejection]

[Date of requesting appeal against examiner's
decision of rejection]

[Date of extinction of right]

Copyright (C); 1998,2003 Japan Patent Office